

**DIGITAL  
IMAGING  
SOLUTIONS**

## **MYVET IMAGING i72W X-Ray Table**

### **The BEST Veterinary Table in the Market**



**Subject:**

**We are the Global Digital X-Ray  
Imaging Component & Solution  
Provider with full lineup of veterinary  
X-Ray Detector products.**




November 8, 2017

# Elevating Table

*Table lowers for easy patient loading*

Our motorized table can go as low as 10" off the ground- Think about all the trouble a large 100 lb dog presents to an office staff when you try to get them on a standard table and how much risk they take with the dogs health if you need to sedate them to take the X-Rays.

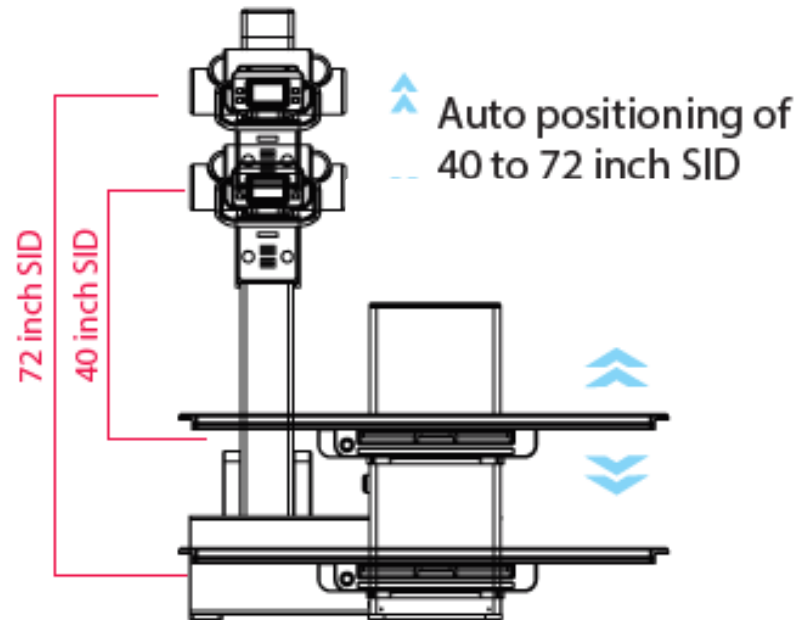

# 40-72 SID

*See more anatomy with the same panel*

Selectable 40-72" SID. The FOV spreads as SID increases. At a 72" SID you will see an added inch of anatomy on all 4 sides due to the decrease in magnification.. Close collimator shutters to reduce x-ray beam spread. Helps reduce scatter radiation and lowers excess exposure. More parallel portion of x-ray beam used that reduces penumbra. A 72" SID avoids the need for an added exposure with more anatomy in a single view and sharper image quality.



40 SID




72 SID

# Image Stitching

*Easily stitches 2 and 3 images*


We have stitching software - What an advantage to be able to take a set of 3 exposures and get the full body of a large dog without the need to move them at all.

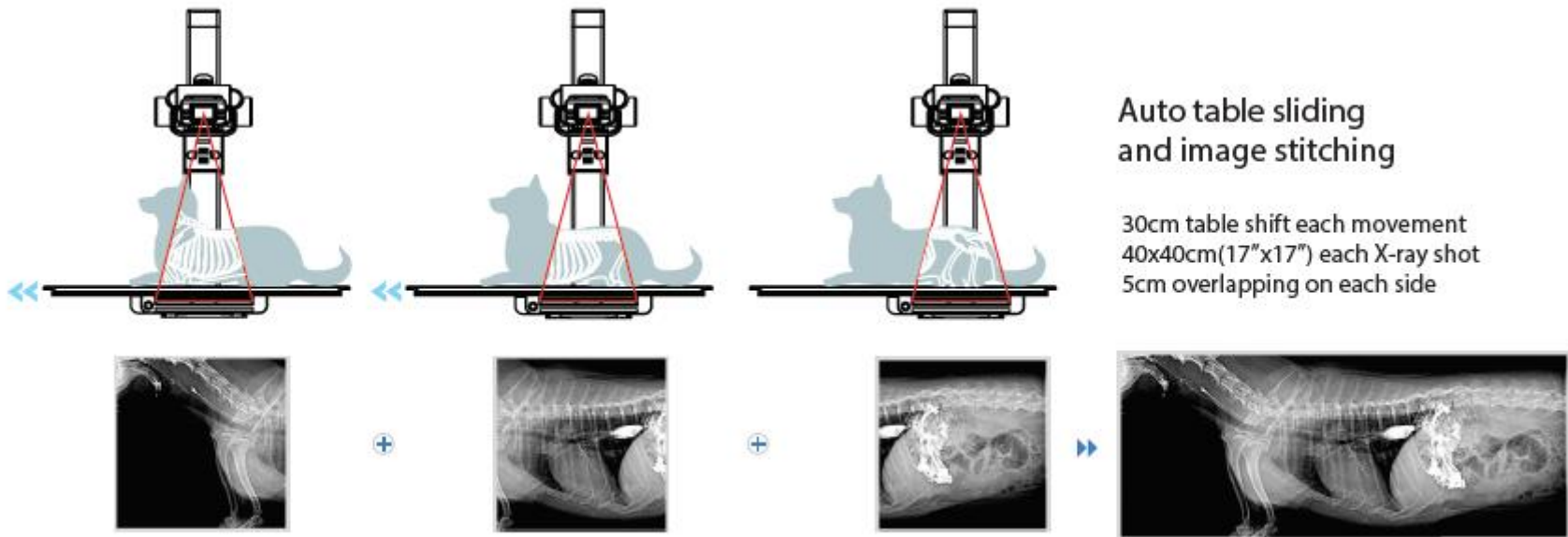

## Weight Bearing

### *Gravity Dependent Images*


We offer a Wall Stand for a 72" SID. Now you can do weight bearing, fluid exams and standing exams reducing pain or discomfort for the animal as you don't need to put them on the table..

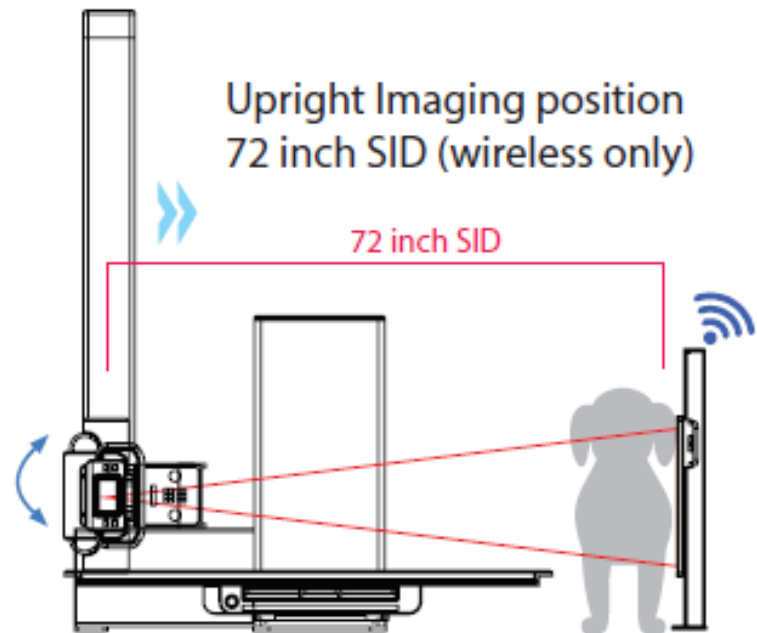

## Quiet Elevation

### *Noise Free Electric Motors for Patient Comfort*


Unique design for Quiet/ Smooth table elevation and table movement- The animal is already in stress and the alien noise of the average table only makes it harder to keep them in position as you need to move them.

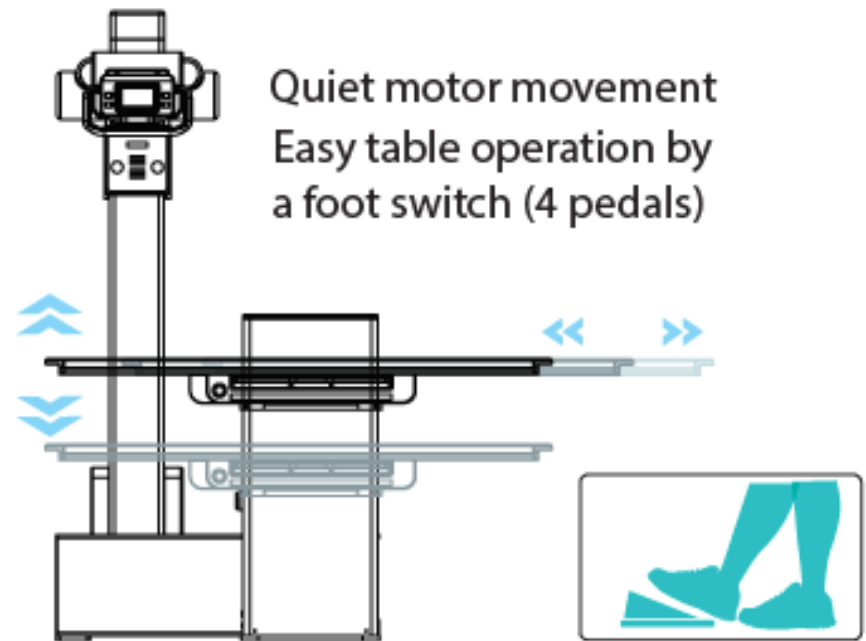

## 40 SID sync

### *Tube and Table move in unison*

Auto tracking of table and tube stand height at 40"- No need to adjust SID as you elevate table which saves time and allows staff to focus on the patient and positioning. This saves time and avoids retakes

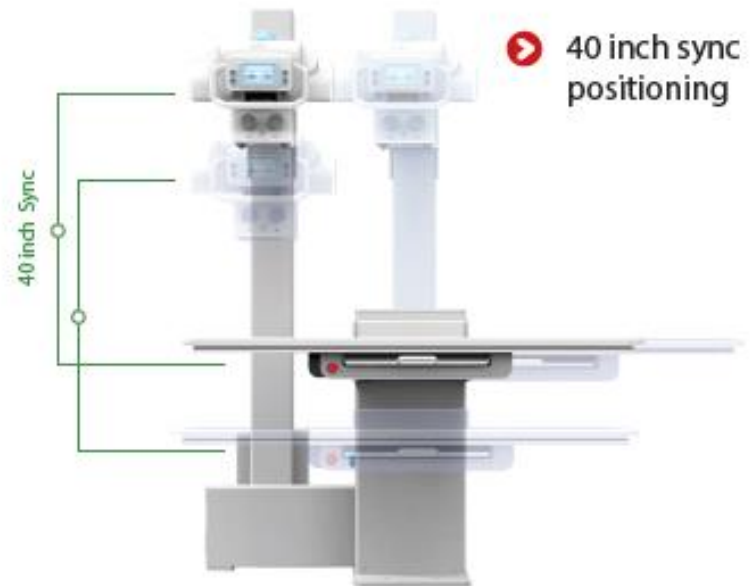


# Bucky Charger

*Wireless FPD that is always ready when needed*

1717 Wireless DR panel with bucky charger- Get the maximum area of coverage, cable free movement from table to wall stand without the need to worry about changing batteries. The 1717WCC is a Cesium panel for highest image resolution and is water resistant with a 600 Lb weight bearing capacity.



 Charging Bucky with a magnetic contact



# Warranty

*Industry Best*

5 years of parts warranty on entire system- No unexpected expenses for the customer which adds to the cost of ownership value plus you get 5 years of Rayence Phone Technical Support with updates



# Support

*One phone call for all support*

Integrated X-Ray System-

RayenceX-Ray+Rayence DR + Rayence Software all integrated- The complete package from one supplier which makes it easy to operate and provide support when needed. One number to call for all service issues.

